

Medtronic


Shiley™ oral/nasal endotracheal tube with intermediate cuff

For paediatric through adult patients.

The Shiley™ oral/nasal endotracheal tube with intermediate cuff can be used for paediatric through adult patients.

- 5.0-10.0 mm size options feature an intermediate high-volume, low-pressure cuff
- 3.0-4.5 mm size options feature an intermediate low-profile cuff to meet the unique needs of paediatric patients

This product is made with non-DEHP PVC material, which softens at body temperature and molds to the airway.

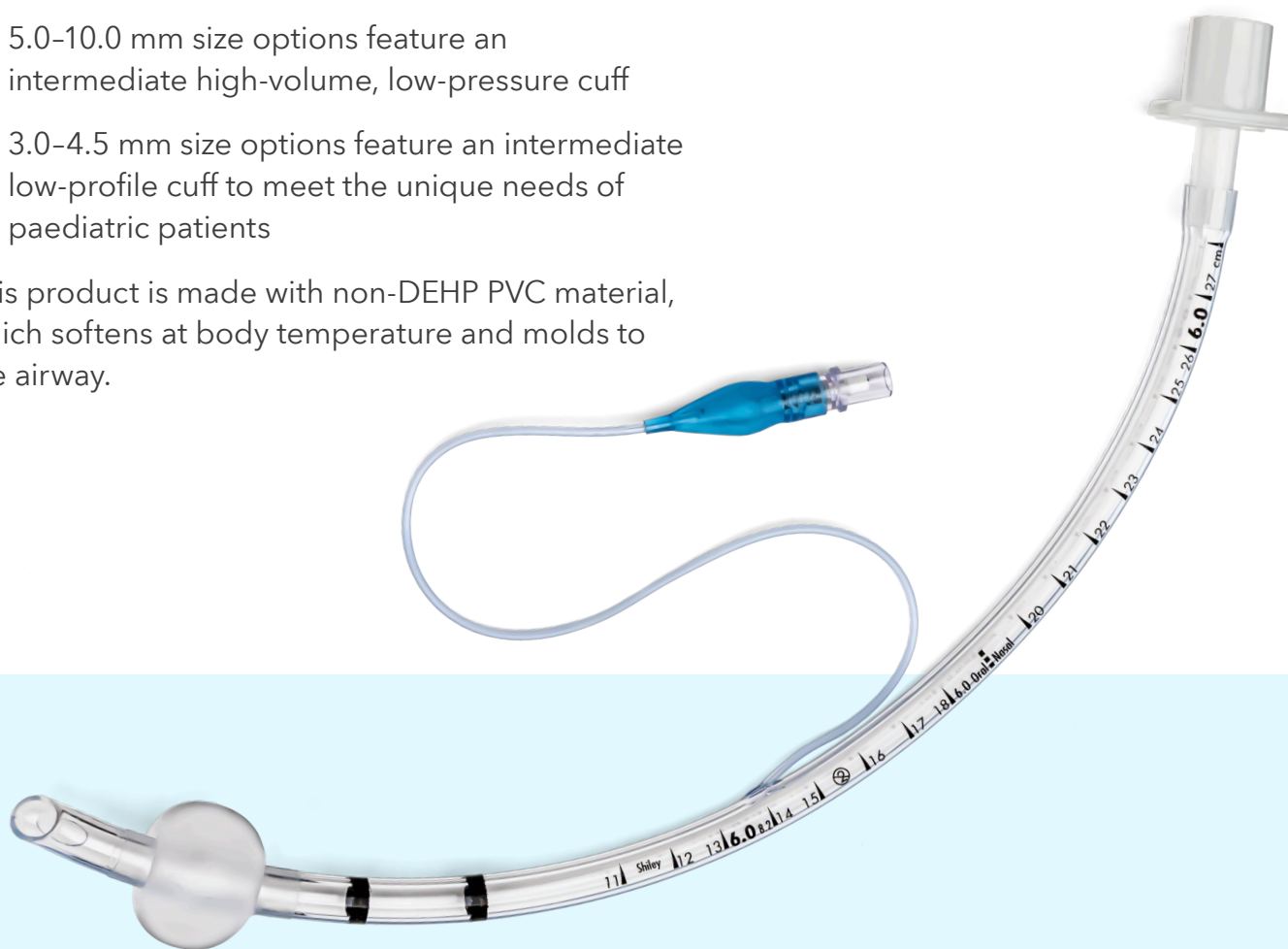

Features and benefits

Sizes 5.0-10.0

A standard 15 mm connector
For connection to respiratory and anaesthesia equipment

Depth marks in centimetre increments
Printed along a portion of the length of the device

Magill curve
For easier tube insertion

Embedded radiopaque filament within its wall
Present down the length of the tube and visible on x-ray

Glottic print marks – full circumference
For tube visualisation during intubation and tube placement when the tip is out of view and beyond the vocal cords†

Intermediate high-volume, low pressure cuff

Hooded tip
To reduce tip sharpness

Beveled end with Murphy eye
To aid in visibility and insertion through the vocal cords and into the trachea

Sizes 3.0-4.5

Glottic print marks – 180°

Intermediate low profile cuff
Elongated, low-profile cuff used in paediatrics

Straight tip with a Murphy eye

Ordering Information

| CFN | I.D. (mm) | O.D. (mm) | Length (mm) | Cuff Ø (mm) |
|-------|-----------|-----------|-------------|-------------|
| 87430 | 3.0 | 4.4 | 160 | 7.8 |
| 87435 | 3.5 | 5.0 | 180 | 9.6 |
| 87440 | 4.0 | 5.7 | 200 | 9.6 |
| 87445 | 4.5 | 6.3 | 220 | 13.0 |
| 87450 | 5.0 | 6.9 | 240 | 20.9 |
| 87455 | 5.5 | 7.5 | 270 | 21.8 |
| 87460 | 6.0 | 8.2 | 280 | 23.9 |
| 87465 | 6.5 | 8.9 | 290 | 24.7 |
| 87470 | 7.0 | 9.5 | 300 | 28.2 |
| 87475 | 7.5 | 10.2 | 310 | 28.7 |
| 87480 | 8.0 | 10.8 | 320 | 32.3 |
| 87485 | 8.5 | 11.4 | 330 | 32.9 |
| 87490 | 9.0 | 12.1 | 335 | 34.0 |
| 87495 | 9.5 | 12.8 | 335 | 35.1 |
| 87410 | 10.0 | 13.5 | 335 | 35.2 |

Carton quantity: 10. All tubes are delivered sterile packed, for single use only.

†As indicated by ISO 5361:2016.

Medtronic

Medtronic Australasia Pty Ltd
2 Alma Road, Macquarie Park,
NSW 2113 Australia
Tel: +61 2 9857 9000
Fax: +61 2 9889 5167
Toll Free: 1800 668 670

Medtronic New Zealand Ltd
Level 3 - Building 5, Central Park Corporate Centre
666 Great South Road
Penrose, Auckland 1051
New Zealand
Fax: +64 9 918 3742
Toll Free: 0800 377 807

medtronic.com.au
medtronic.co.nz

Please refer to the product insert or Instructions for Use for a complete list of indications, contraindications, warnings, precautions, and other important medical information. Always follow the Instructions for Use.

© Medtronic 2022 All Rights Reserved
RI 007-11-22 ANZ. #11958-112022