


A transparent incise film which acts as a bacterial barrier

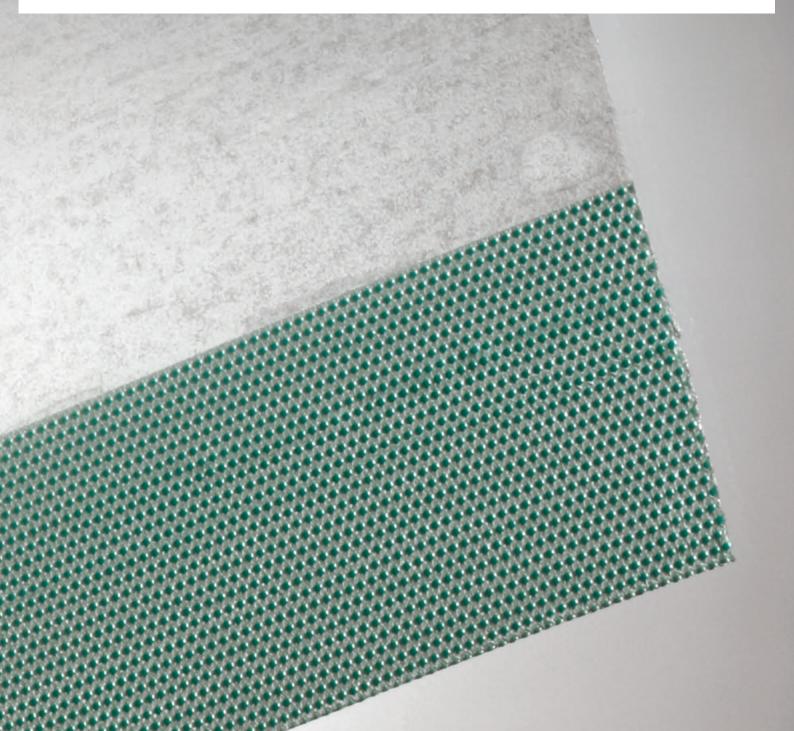

OPSITE^{*} Incise is a protective film that is a barrier to bacteria

OPSITE Incise is a transparent, adhesive polyurethane film. It allows the skin to breathe and therefore prevents the build up of moisture under the drape. As a result OPSITE Incise adheres throughout long surgical operations to the surrounding skin and most importantly to the wound edge. This means that it also helps to prevent the lateral migration of bacteria which can lead to infection.

OPSITE Incise holds theatre drapes securely in place. They also provide a sterile field around the incision site.

Advantages

Adheres right up to wound edge

The drape will provide protection from laterally migrating bacteria because it stays in place even throughout major operations involving prolonged retraction of wound edges.

Transparent minimal glare surface.

Thin yet strong

OPSITE Incise is strong yet can be cleanly cut with a scalpel

Conformable and elastic

OPSITE Incise is very conformable on awkward operation sites, yet its elastic qualities make it suitable for use where heavy manipulative surgery is involved, for example, in a total hip replacement.

Maintains a sterile operation site

OPSITE Incise Drape provides a sterile working field.

Indications

Suitable for all types of surgery:

- 1. Orthopaedic
- 2. Neuro surgery
- 3. Major abdominal
- 4. Plastic
- 5. Open Heart and Thoracic
- 6. Ophthalmic
- 7. Paediatric

Note on application

The skin should be thoroughly prepared in the normal way. Before applying OPSITE Incise, the skin should be allowed to dry after using cleaning solutions e.g. iodine.

Presentation

Each OPSITE Incise Drape is sterile and individually wrapped in a peel apart pouch.

Sterilisation is by ethylene oxide.

Ordering OPSITE Incise

Code	Product description	Items per unit
4967	14cm x 25cm	Box/20*
4986	28cm x 15cm	Box/10
4987	28cm x 30cm	Box/10
4988	28cm x 45cm	Box/10
4995	42cm x 40cm	Box/10
4989	55cm x 45cm	Box/10
4994	84cm x 56cm	Box/10

*Not available in New Zealand

Please consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.

 Smith & Nephew Pty Ltd

 Australia

 T 61 2 9857 3999

 F 61 2 9857 3990

 www.smith-nephew.com/australia

 Customer Service

 T 13 13 60

 F 1800 671 000

 New Zealand

 T 64 9 820 2840
 F 64 9 820 2841

 www.smith-nephew.com/new-zealand

 Customer Service

 T 0800 807 663
 F 0800 263 222







