


## NEW PRODUCT IN OUR RANGE

---

### ANAL SAC EXCISION KIT • NEW GENERATION

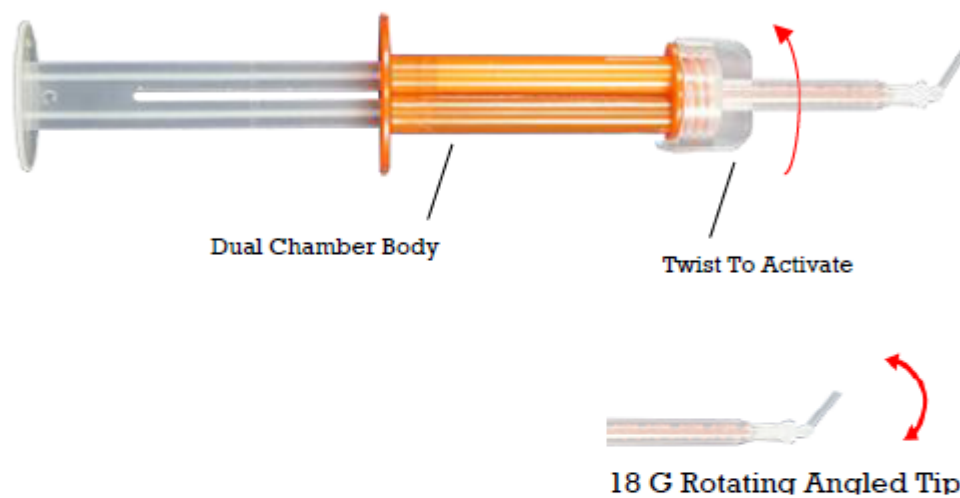

New generation of anal sac excision kit for easier and quicker operation

- Dual chambered cartridge each with a nontoxic inert polymer
- Mixing tip: the 2 components are mixed when pushed in the mixing chamber
- 18 G plastic blunt angled needle for easy insertion into anal sac opening
- The gel solidifies 2-4 minutes upon application
- Capacity of 2.4 ml, sufficient for larger patients
- Ergonomic design for small hands
- For activation: twist the distal part clockwise till a click is heard